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Histological slide viewer, user guide

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Conexión a internet: mínimo 20mb/s

Navegador: Google Chrome, Firefox

Welcome to the histological slide viewer. This brief tutorial will help you navigate the viewer easily. We have aimed to provide a simple and intuitive tool for viewing histological slides via a web browser, so it can be accessed from any computer (windows/mac/linux) with an internet connection. To ensure a comfortable and speedy experience, we recommend using devices that meet these minimum requirements:

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
To start slide visualization follow the steps

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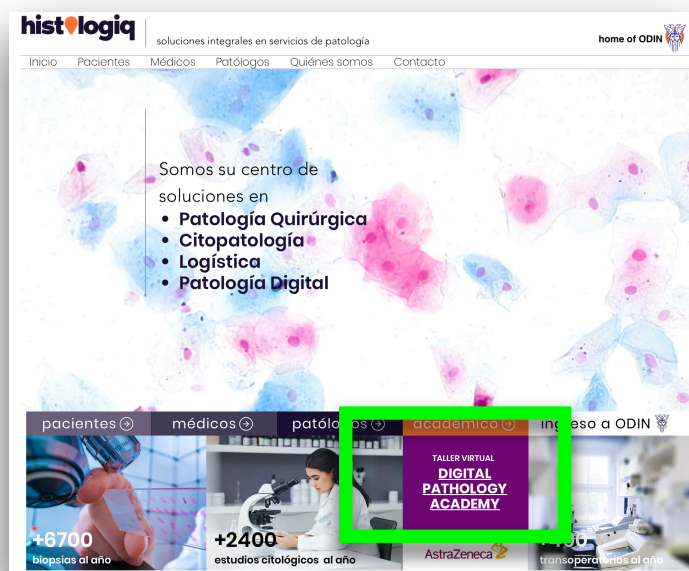Once selected, the system will display the list of slides available for evaluation.

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Once selected the desired slide, the system will ask for username and password (unique for each webinar), which is annotated in the header of the webpage (see pink arrow)

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


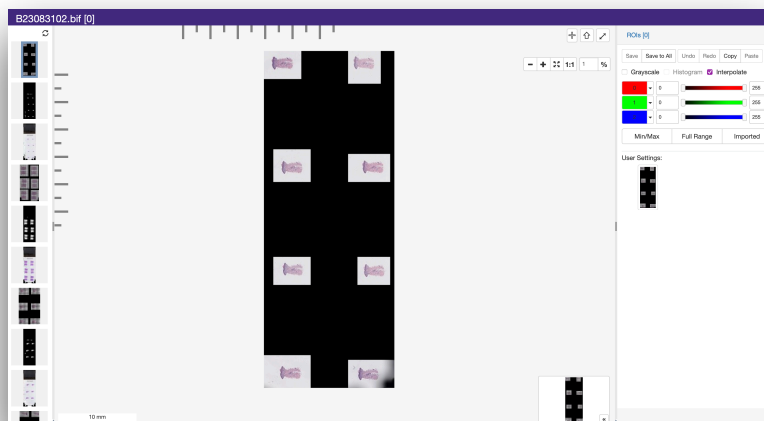
The slide viewer will provide a panoramic image of the slide. For optimization of the browsing experience, only the areas with tissue are represented (areas of the slide without tissue are removed by the program and appear black). The different components of the viewer are specified in the following image

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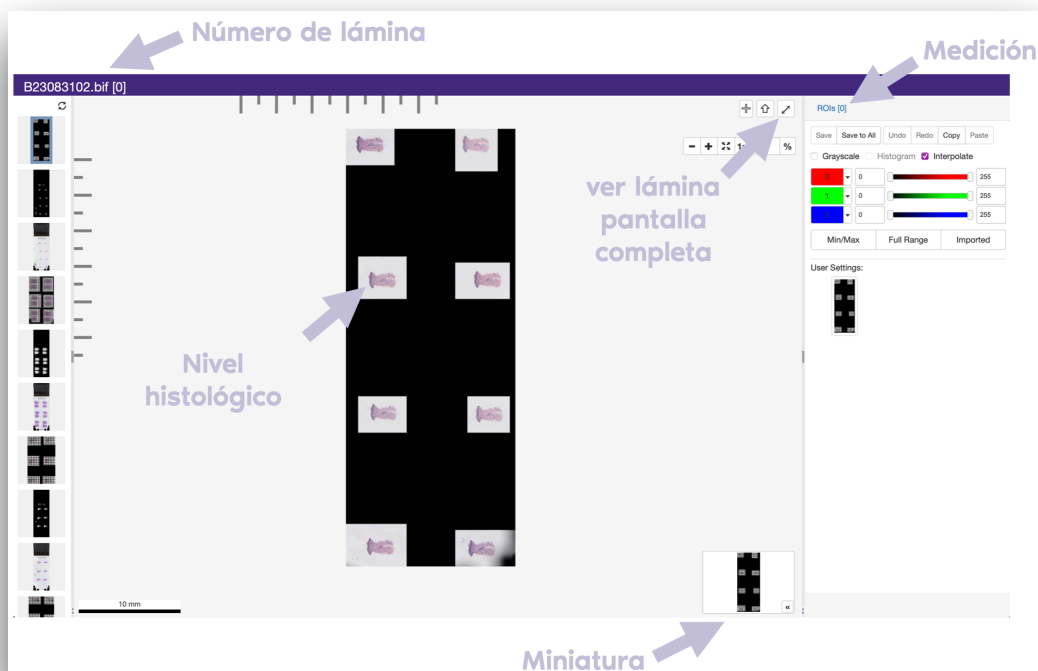To zoom in or out of the image, you can either press the - and + buttons in the upper right corner or use the mouse (with the scroll wheel or a Magic Mouse/trackpad)

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



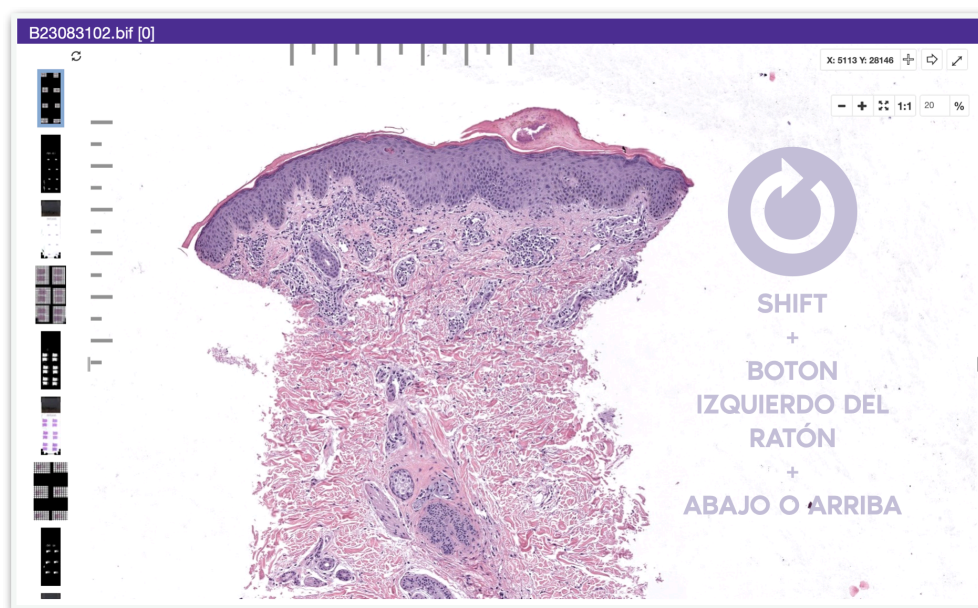
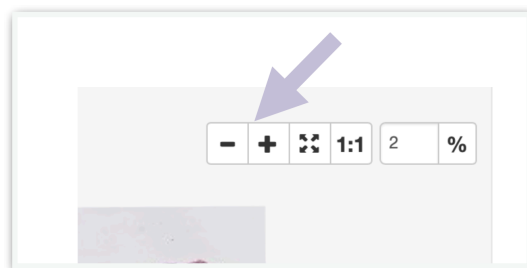
To rotate the slide, simultaneously press the SHIFT key, the left mouse button, and move it up or down as desired

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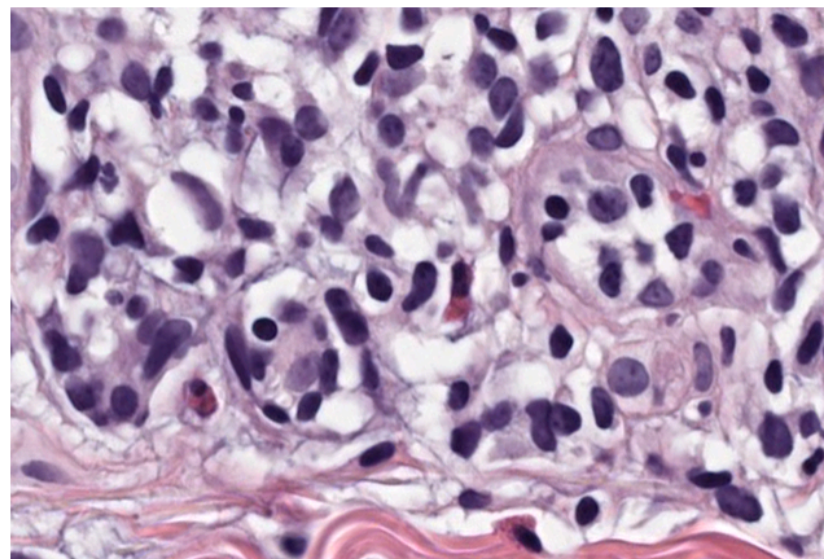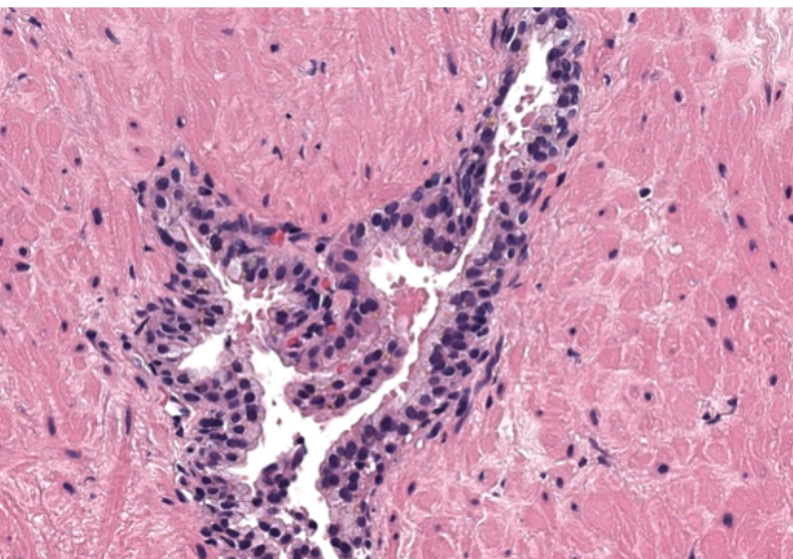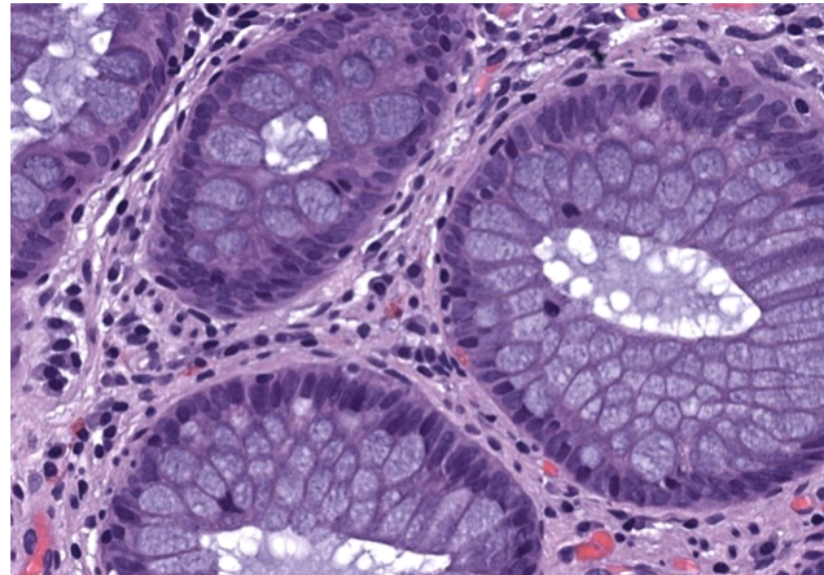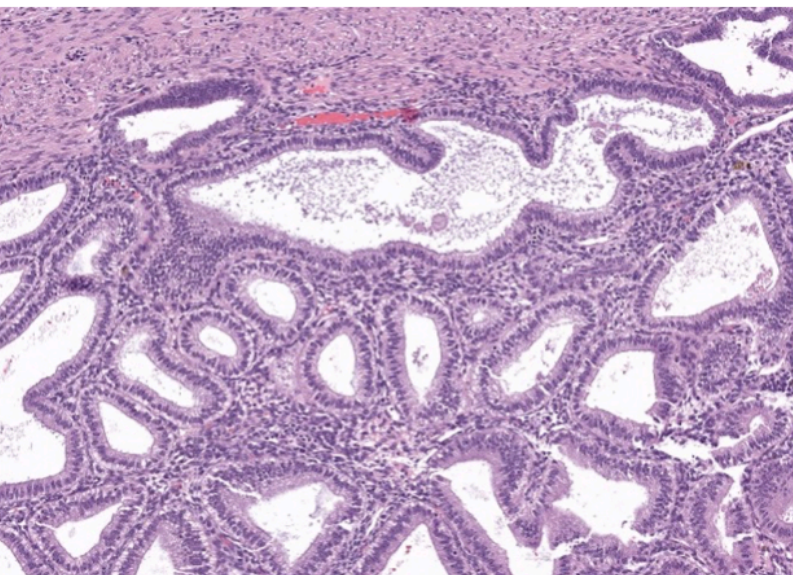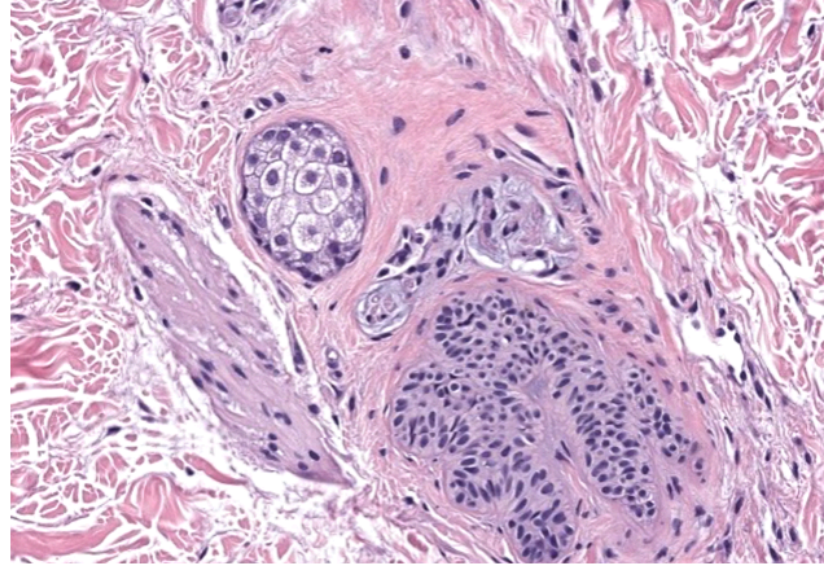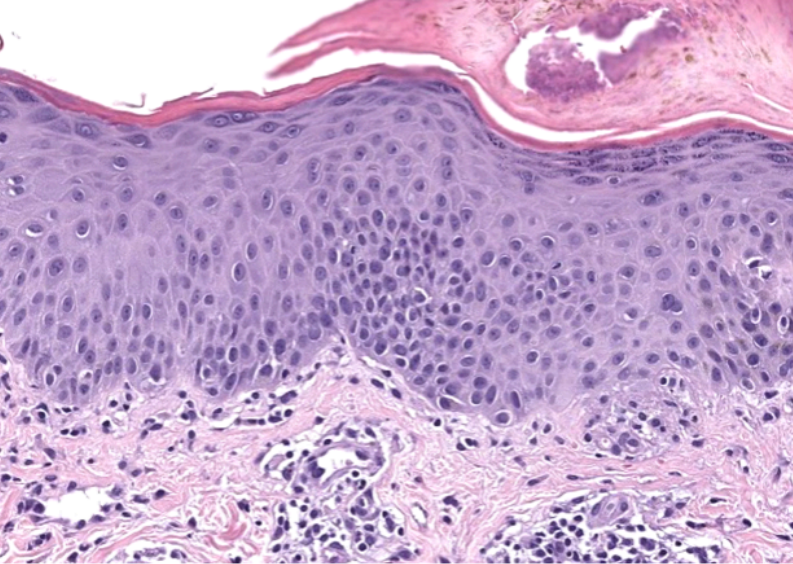We are convinced that digital pathology is the new standard in the practice of surgical pathology, and we want to make it a pleasant and user-friendly experience. Please do not hesitate to contact us if you need guidance on using this tool.

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